

ZytoLight

SPEC EGFR/CEN 7 Dual Color Probe

REF Z-2033-200 ∇_{Σ} 20 (0.2 ml)

REF Z-2033-50 ∇_{Σ} 5 (0.05 ml)


For the detection of the human EGFR gene and alpha-satellites of chromosome 7 by fluorescence *in situ* hybridization (FISH)



IVD

In vitro diagnostic medical device

according to EU directive 98/79/EC

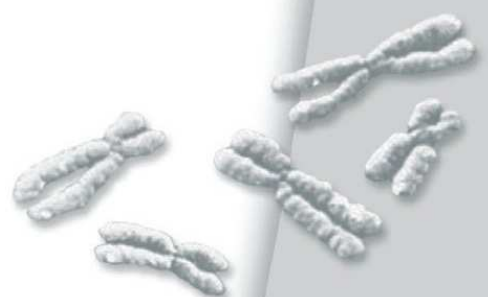

Fluorescence-labeled polynucleotide probe for the detection of the human EGFR gene and alpha-satellites of chromosome 7 centromeres, ready to use

Product Description

- Content:** ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15) in hybridization buffer. The probe contains green-labeled polynucleotides (ZyGreen: excitation at 503 nm and emission at 528 nm, similar to FITC), which target the EGFR gene, and orange-labeled polynucleotides (ZyOrange: excitation at 547 nm and emission at 572 nm, similar to rhodamine), which target alpha-satellite-sequences of the centromere of chromosome 7.
- Product:** Z-2033-200: 0.2 ml (20 reactions of 10 µl each)
Z-2033-50: 0.05 ml (5 reactions of 10 µl each)
- Specificity:** The ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15) is designed to be used for the detection of the human EGFR gene as well as chromosome 7 alpha-satellites in formalin-fixed, paraffin-embedded tissue or cells by fluorescence *in situ* hybridization (FISH).
- Storage/Stability:** The ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15) must be stored at 2...8°C protected from light and is stable through the expiry date printed on the label.
- Use:** This product is designed for *in vitro* diagnostic use (according to EU directive 98/79/EC). Interpretation of results must be made within the context of the patient's clinical history with respect to further clinical and pathologic data of the patient by a qualified pathologist!

Safety Precautions:

Read the operating instructions prior to use!

Do not use the reagents after the expiry date has been reached!

This product contains substances (in low concentrations and volumes) that are harmful to health. Avoid any direct contact with the reagents. Take appropriate protective measures (use disposable gloves, protective glasses, and lab garments)!

If reagents come into contact with skin, rinse skin immediately with copious quantities of water!

A material safety data sheet is available on request for the professional user!

Principle of the Method

The presence of certain nucleic acid sequences in cells or tissue can be detected by *in situ* hybridization using labeled DNA probes. The hybridization results in duplex formation of sequences present in the test object with the labeled DNA probe.

Duplex formation (with sequences of the EGFR gene and chromosome 7 alpha-satellites in the test material) is directly detected by using the tags of fluorescence-labeled polynucleotides.

Instructions

Pretreatment (dewaxing, proteolysis, post-fixation) should be carried out according to the needs of the user.

Denaturation and hybridization of probe:

- 1.** Pipette 10 μ l ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15) each onto individual samples

A gentle warming of the probe, as well as using a pipette tip which has been cut off to increase the size of the opening, can make the pipetting process easier. Avoid long exposure of the probe to light.

- 2.** Avoiding trapped bubbles, cover the samples with a coverslip (22 mm x 22 mm). Seal the coverslip, e.g. with a layer of hot glue from an adhesive pistol or with rubber cement

- 3.** Denature the slides at 75°C ($\pm 2^\circ$ C) for 10 min, e.g. on a hot plate

Depending upon the age of the sample and variations in the fixation stage, it may be necessary to optimize the denaturing temperature (73°C-77°C).

- 4.** Transfer the slide to a humidity chamber and hybridize overnight at 37°C (e.g. in a hybridization oven)

It is essential that the tissue/cell samples do not dry out during the hybridization step.

Further processing, such as washing and counter-staining, can be completed according to the user's needs. For a particularly user-friendly performance, we recommend the use of a ZytoLight FISH system by ZytoVision. These systems were also used for the confirmation of appropriateness of the ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15).

Results

With the use of appropriate filter sets, the hybridization signals of the labeled EGFR gene appear green; the hybridization signals of the labeled alpha-satellite-sequences of the centromere of chromosome 7 appear orange. In interphases of normal cells or cells without aberrations of chromosome 7, two EGFR signals and two chromosome 7 signals appear. In cells with a gene amplification an increased number of gene specific signals or signal clusters are visible.

The polynucleotides contained in the ZytoLight SPEC EGFR/CEN 7 Dual Color Probe (PL15) which recognize the alpha-satellite-sequences of the centromere of chromosome 7 function in themselves as an internal control that a successful hybridization has occurred, as well as proving the integrity of the cellular DNA.

In order to judge the specificity of the signals, every hybridization should be accompanied by controls. We recommend using at least one control sample in which the chromosome 7 and EGFR gene region copy number is known.

Care should be taken not to evaluate overlapping cells, in order to avoid false results, e.g. an amplification of genes. Due to decondensed chromatin, single FISH signals can appear as small signal clusters. Thus, two or three signals of the same size, separated by a distance equal to or less than the diameter of one signal, should be counted as one signal.

Our experts are available to answer your questions.

Literature

- Alexandrov IA, et al. (1988) *Chromosoma* **96**: 443-53.
- Cappuzzo F, et al. (2005) *Natl Cancer Inst* **97**: 643-55.
- Kievits T, et al. (1990) *Cytogenet Cell Genet* **53**: 134-6.
- Kondo I, Shimizu N (1983) *Cytogenet Cell Genet* **35**: 9-14.
- Libermann TA, et al. (1985) *Cell Sci Suppl* **3**: 161-72.
- Merlino GT, et. (1985) *J Clin Invest* **75**: 1077-9.
- Sassen A, et al. (2008) *Breast Cancer Res* **10**: R2.
- Tovey SM, et al. (2004) *Breast Res* **6**: 246-51.
- Waye JS, Willard HF (1986) *Nucleic Acids Res* **14**: 6915-27.
- Wilkinson DG: In Situ Hybridization, A Practical Approach, *Oxford University Press* (1992)
ISBN 0 19 963327 4.

As of: June 10, 2013 (5.1)

Trademarks:

ZytoVision® and ZytoLight® are trademarks of ZytoVision GmbH.



ZytoVision GmbH · Fischkai 1


D - 27572 Bremerhaven · Germany

Phone: +49 (0) 471/4832 - 300

Fax: +49 (0) 471/4832 - 509

www.zytovision.com

info@zytovision.com

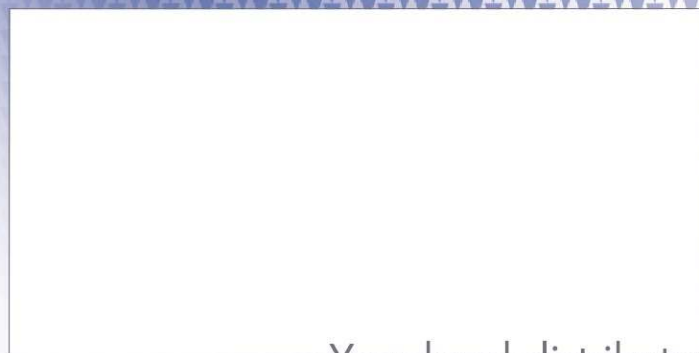

Your local distributor